

sagemax.com

sagemax®



Sagemax experience

---

NexxZr T Multi monolitna  
kruna vs slojevanje keramike  
(cut-back) by Ivan Unfeter

NexxZr T Multi monolitna kruna vs slojevanje keramike (cut-back)  
by Ivan Unfeter

# Sagemax experience


## Šema karakterizacije

## Šema slojevanja keramike

**Materijal**  
IPS Ivocolor, Co. Ivoclar


### Step 1

**Boja/Mikstura**  
SD1 + E02 Cream +  
Glaze Fluo  
**Odnos**  
2/1/2

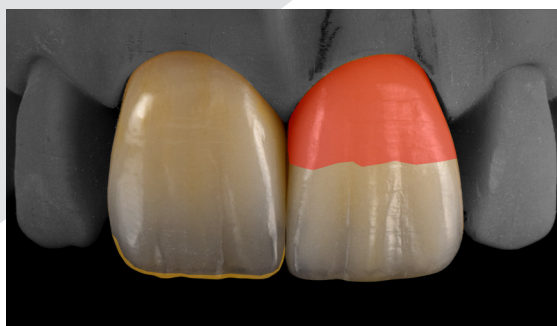


### Step 2

**Boja/Mikstura**  
E15 Ocean + E18 Black  
+ Glaze Fluo  
**Odnos**  
1/1/1

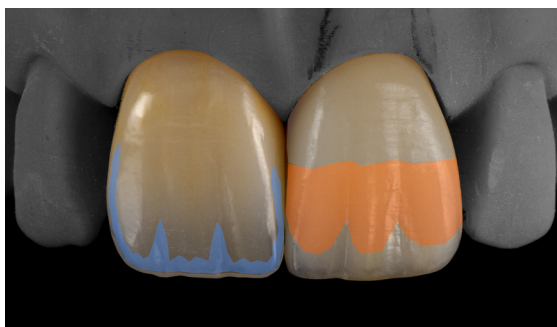


### Step 3

**Boja/Mikstura**  
E02 Cream + SD1 +  
Glaze Fluo  
**Odnos**  
1/1/1

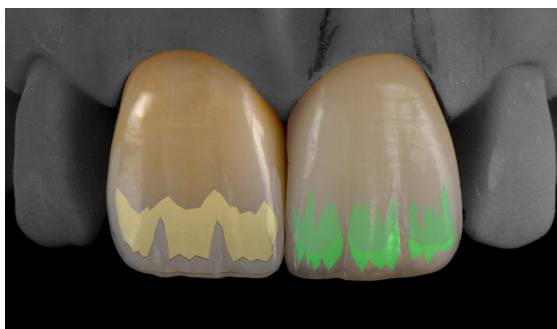

### Step 4

**Boja/Mikstura**  
SD1 + Glaze Fluo  
**Odnos**  
1/1

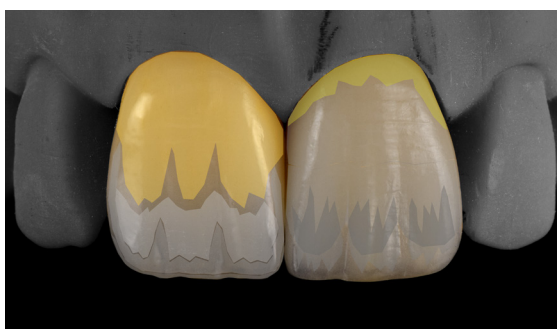

**Materijal**  
IPS e.max Ceram, Co. Ivoclar

### Step 1

**Boja/Mikstura**  
Power Dentin A1  
**Odnos**  
Čisto

### Step 2

**Boja/Mikstura**  
Dentin A1  
**Odnos**  
Čisto

### Step 3

**Boja/Mikstura**  
Dentin A1 +  
Mamelon light  
**Odnos**  
1/1

🔥 Pečenje dentina

### Step 4

**Boja/Mikstura**  
Honey  
**Odnos**  
Čisto

## Step 5

**Boja/Mikstura**  
E01 White + Glaze Fluo  
**Odnos**  
1/1

Stain pečenje 🔥




## Step 5

**Boja/Mikstura**  
OE 2  
**Odnos**  
Čisto

## Step 6


**Boja/Mikstura**  
Glaze Fluo  
**Odnos**  
Čisto

Pečenje glazure 🔥

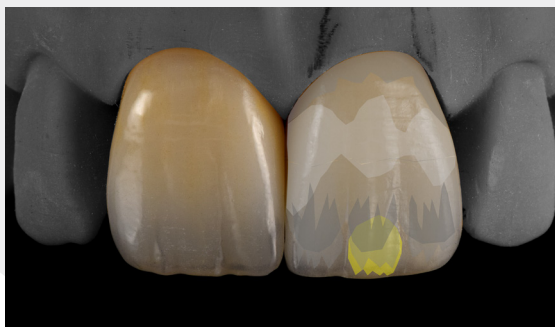


## Step 6

**Boja/Mikstura**  
Citrine  
**Odnos**  
Čisto

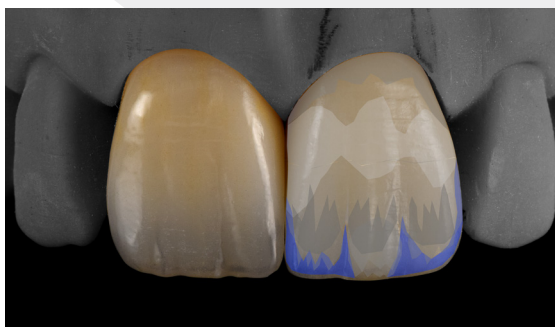

## Step 7


**Boja/Mikstura**  
OE 1  
**Ratio**  
Čisto

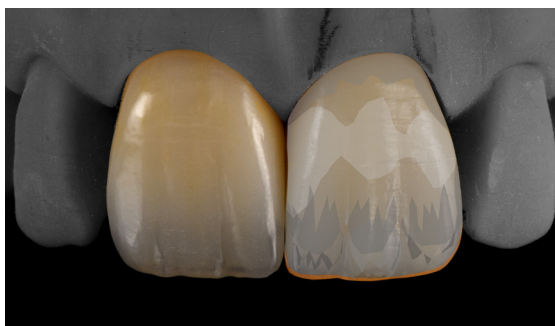

## Step 8

**Boja/Mikstura**  
Dentin A1  
**Odnos**  
Čisto

🔥 Enamel pečenje

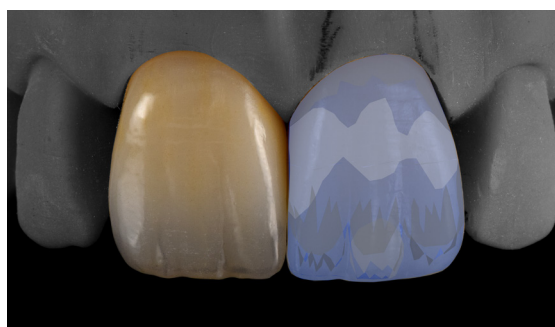


## Step 9

**Boja/Mikstura**  
Transpa Incisal 1  
**Odnos**  
Čisto

🔥 Enamel/Korektivno  
pečenje

## Finalni rezultat

---

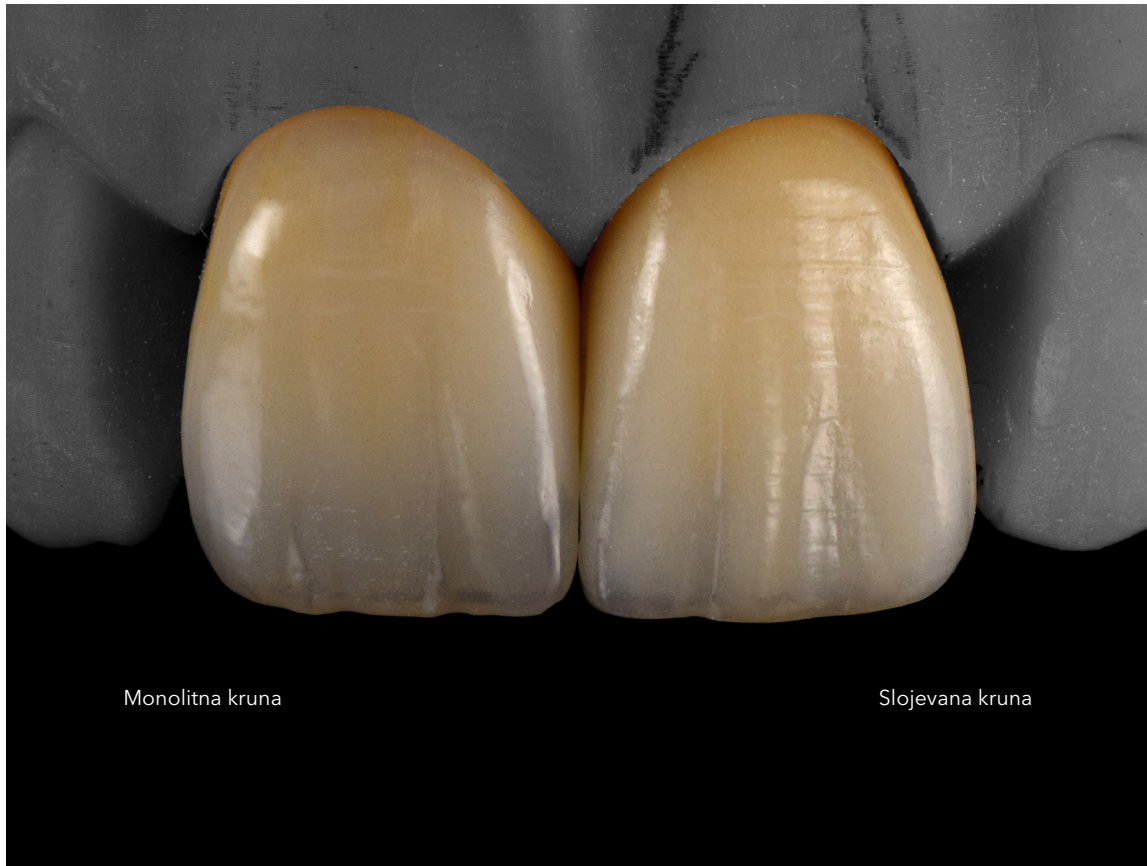



Monolitna krana

Slojevana krana



### Dental designer

**Ivan Unfeter**  
Dental Technician, Serbia

 dent.artist.unfeter  
 dent.artist.unfeter

---

Contact us  
[info@sagemax.com](mailto:info@sagemax.com)

[sagemax.com](http://sagemax.com)

Follow us  
[@sagemax.official](https://www.instagram.com/sagemax.official)

